


Unrestricted Access and Stability

Mutilok™
Hand Surgery
Table System

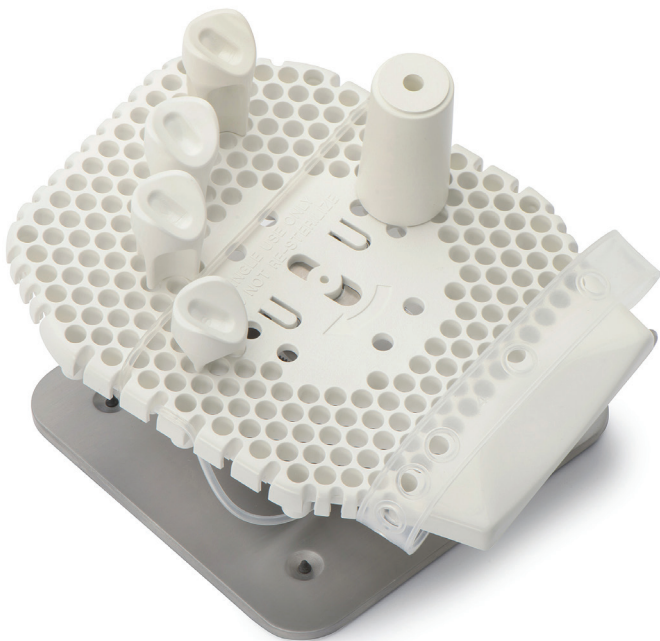

Multilok™ Hand Surgery Table System

A Surgical Hand Positioner for Simplified Procedures

The **Multilok™ Hand Surgery Table System** represents a major advance in operating systems for general and micro-surgical procedures of the wrist, hand and fingers. The surgeon is able to select from a range of finger posts and silicone straps, a combination of components that will give the most appropriate hand position.

Features

- Range of accessory components offers the surgeon multiple setup options
- Pre-sterilized disposable accessory components for surgical convenience
- Disposable table top and accessories
- Unrestricted access to operative site
- Infinite number of positions are possible
- Simple adjustments made by releasing the unique cam lock lever
- Stainless steel hand rest for improved surgeon comfort and stability
- Pre-sterilized table top and accessories constructed of high quality plastic for single use



Multilok Modular System

The Multilok system consists of four major modular components:

1. A resterilizable stainless steel base with built in table locking mechanism*
2. A stainless steel hand rest to support and steady the surgeon's hand at the appropriate height
3. A pre-sterilized, polished table top, constructed of high quality plastic for single use
4. A range of finger supports, thumb posts, silicone straps and silicone retractors

* Attaches to any general surgical table.

Ordering Information

REF	
0302-000	Multilok Starter Set (1 each of 0302-001 through 0302-009-10)
0302-001	Multilok Table Base
0302-002	Multilok Table Top (10/cs)
0302-003	Multilok Finger Post, Small (10/cs)
0302-004	Multilok Finger Post, Medium (10/cs)
0302-005	Multilok Finger Post, Large (10/cs)
0302-006	Multilok Finger Post, Thumb (10/cs)
0302-007	Multilok Wrist Strap (10/cs)
0302-008	Multilok Retractor (4/cs)
0302-009-10	Multilok Silicone Ties (10/cs)
0302-010	Multilok Hand Rest
0302-012	Multilok Blunt Silicone Retractor (4/cs)



30031 Ahern Avenue
Union City, CA 94587-1234 USA
Telephone: 510-429-1500
Toll Free: 800-777-4674
Fax: 510-429-8500
Outside USA: +1-510-429-1500
mizuho.com
newhipnews.com



Corin Group PLC
The Corinium Center
Cirencester, GL71YJ
United Kingdom
Tel: +44 (0) 1285 659 866



Note: Mizuho OSI is constantly improving its products. All specifications are subject to change without notice. Multilok™ is a trademark of Corin Medical Ltd. The Multilok™ Hand Operating Table is distributed in the USA on behalf of Corin Medical Ltd. • Mizuho OSI is a Delaware Corporation.

2016 ©Mizuho OSI

REF NW0890 REV A